Veiligheid

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Safety

To ensure long-lasting and optimal use of the Information Speaker BT, we recommend reading this manual carefully before using the speaker. This will help prevent damage caused by incorrect use. Pay particular attention to the safety instructions. If you transfer the Information Speaker BT to someone else, please also pass on this manual. Keep the manual for future reference.

Safety Instructions

Do not perform any disassembly or assembly work on the device. Do not rinse the device under running water and do not submerge it in water. Do not expose the speaker to extreme heat (sun, heaters). Stop using the device immediately if smoke, odor, leakage, or deformation occurs. Avoid short-circuiting the USB port. Store the speaker in a dry place, preferably in its original packaging, to prevent moisture damage. The speaker is not drop- or shock-resistant. During prolonged use at high volume, the device may become warm. This is normal. The speaker will cool down to its normal temperature once the volume is reduced.

Introduction

The Information Speaker BT is a sound exciter with a carbon vibration plate that is placed directly on the skin. The vibrations turn the body into a resonant chamber, allowing frequencies to be not only heard but also felt deep within the body. The speaker combines tactile vibrations, audible tones, and electromagnetic light frequencies for optimal information transfer.

DAMP


a combination of melanin, gold, and carbon - is integrated into the speaker and supports the transmission of harmonic frequencies. This makes the Information Speaker fundamentally different from standard speakers.

The Information Speaker is designed for use on the body and brings together audible and tangible frequencies in a compact, user-friendly form for optimal information transfer.

Package Contents

Information Speaker BT
Power bank with charging cable
USB-C power cable
Cotton pouch
Manual

Overview

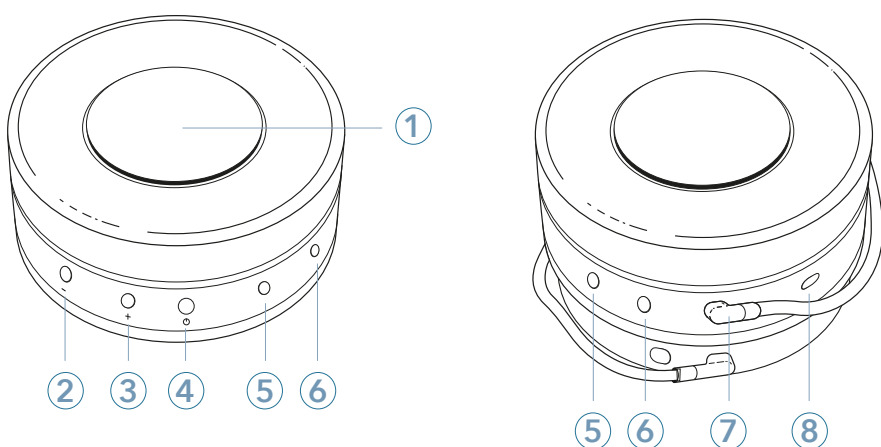1. Carbon vibration plate (to place on the skin)
2. Volume down (VOL-)
3. Volume up (VOL+)
4. Power button (ON/OFF)
5. Red LED - Bluetooth active
6. Green LED - power on
7. USB-C power input
8. Mini-jack (3.5 mm) audio input

Powering on and off: press and hold ON/OFF for several seconds.

Using the Speaker Without Bluetooth

When the audio cable is connected, Bluetooth automatically switches off. The power bank can be used optionally. When the power bank is connected, vibrations and audio output are stronger.

1. Connect your device to the Information Speaker BT using the 3.5 mm audio cable.
2. Optional: Connect the power bank using the USB-C power cable.

Bluetooth Connection

The Information Speaker is set to Single Mode by default (1 speaker).

1. Connect the Information Speaker BT to the power bank using the supplied USB-C power cable. The green LED indicates that the speaker is powered on.
2. Switch on both devices. The red LED indicates that Bluetooth is active.
3. On your phone, tablet, or laptop, search for "info med speaker" and pair the device.
4. If pairing does not occur automatically, briefly press VOL+ and VOL- at the same time and accept the pairing request.

Tips

- If <<info med speaker>> does not appear in the list, briefly switch Bluetooth off and on again on your playback device. Make sure your device actively scans for available connections.
- Depending on your playback device, it may take a moment for the speaker to appear – a little patience works wonders.

Pairing two speakers in sync (*Twin-modus*)

Step 1: Pair the first speaker

1. Press VOL+ and VOL- simultaneously on Speaker 1.
2. Accept the pairing request on your phone or tablet.

Step 2: Set both speakers to Twin Mode

- Speaker 1:**
1. Press ON/OFF three times.
 2. Within 10 seconds, press VOL+.
- Speaker 2:**
1. Press ON/OFF three times (LED flashes).
 2. Within 10 seconds, press VOL+.
- ✓ Both speakers are now in Twin Mode.

Step 3: Is Speaker 2 not pairing automatically?

1. Turn both speakers off (hold ON/OFF).
2. Turn Speaker 1 on.
3. Wait for 2 seconds.
4. Turn Speaker 2 on (it will pair automatically).

Done! Your speakers are now paired in sync and will play the same audio.

Tips

- Not working right away? Repeat Step 3.
- Always wait 2 seconds between turning on Speaker 1 and Speaker 2.
- Both speakers must remain in Twin Mode for synchronized playback.

Resetting

If the speaker does not respond or fails to establish a connection, switch it off and on again manually. If the issue persists, toggle Bluetooth off and on again on your device and retry pairing. You can also reset the speaker to the correct mode:

- Single Mode: when using 1 speaker
- Twin Mode: when using 2 or more speakers in sync

Switching Between Single Mode and Twin Mode

- To switch to Single Mode: press ON/OFF three times, then press VOL- within 10 seconds.
- To switch to Twin Mode: press ON/OFF three times, then press VOL+ within 10 seconds.

Volume Advice

You can adjust the volume using the speaker's controls or through your connected device. During playback, low and high tones may alternate; some tones may be barely audible while others are more prominent. This is normal and part of the frequency composition used by the Information Speaker BT. The ideal volume is one where vibrations are clearly felt while the sound remains softly audible and non-irritating. Excessively high volume may interfere with optimal information transfer. Since everybody is different, the optimal volume level will vary from person to person—and from animal to animal.

Power Consumption

After switching off the speaker, a small amount of power is still drawn briefly until the power bank shuts down automatically. When the speaker is powered on, it sends a small keep-alive signal to the power bank to maintain the power supply.

Audio Indicators

- Startup tone - The speaker is turning on
- Two high beeps - Single Mode is active
- Four high beeps - Searching for second speaker (Twin Mode)
- Four low beeps - Successfully connected to the second speaker
- Short beep - Connected to your device

Using the Information Speaker BT

Place the Information Speaker BT on the chest, positioning the vibration plate at the thymus, approximately 2 cm below the notch between the collarbones. The speaker may also be used locally on other areas of the body, such as the lower back or on specific points of discomfort. For the best results, place the vibration plate directly on bare skin. The speaker is not drop- or shock-resistant.

Use with animals

Place the speaker on the tuber sacrale or directly on the area of concern. Ensure that the speaker is positioned securely and cannot fall. For horses, a thin cotton blanket with a small pocket can be used to hold the speaker in place. For calm dogs and cats, the speaker may be placed on or against the body, either directly on the fur or inside a thin cotton pouch. For small animals, the speaker may be positioned next to the body.

For more detailed instructions and user tips, visit:
<https://informationmedicine.org/en/manuals/>

Cleaning

Clean the device with a soft, dry cloth. If needed, a small amount of isopropyl alcohol may be used.

Do not use aggressive cleaning agents, water, soap, or submerge the device. Do not use metal objects to clean the device. Avoid moisture entering any openings or connectors. Damage caused by improper cleaning is not covered by the warranty.

Technical Specifications

The Information Speaker BT is developed and produced in the Netherlands.

- | | |
|---|---|
| ■ Full-range exciter with neodymium motor | ■ Audio input: 3.5 mm jack |
| ■ Coil: 2.54 cm | ■ Power supply: USB-C |
| ■ Output Power: 10 W RMS | ■ Bluetooth: yes, Single & Twin Mode |
| ■ Impedance: 4 Ohm | ■ Output: Mono |
| ■ Dimensions: Ø 76 mm, height 30 mm | ■ Internal antenna: 2.4 GHz, 5 dBi |
| ■ Weight: 182 grams | ■ Built-in amplifier (powered via power bank) |
| ■ Resonance frequency: 306 Hz | |

The supplied power bank can be magnetically attached to the bottom of the speaker. The power bank includes its own quick-start guide. Follow this separate manual for safe operation. **Warning:** Do not use power supplies with an output voltage above 5 volts (fast chargers). Higher voltage may damage the charging cable and create a fire hazard.

Warranty

You have a 24-month warranty on the Information Speaker BT and a 12-month warranty on the audio cable. The product has been thoroughly tested before shipment. The invoice serves as proof of warranty. The warranty does not cover damage caused by drops or impacts, moisture damage, misuse, or unauthorized repairs.

For full warranty conditions, visit <https://informationmedicine.org/en/general-terms-conditions/>.

Recycling and Disposal

The packaging materials are recyclable. Dispose of packaging according to local regulations. This product may not be discarded with household waste within the European Union. Dispose of electronic devices in accordance with regional electronic waste regulations.

Copyright

This document is protected by copyright. Reproduction or reprinting, even in part, and the use of images – even in modified form – is only permitted with written permission from Information Medicine BV.

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